
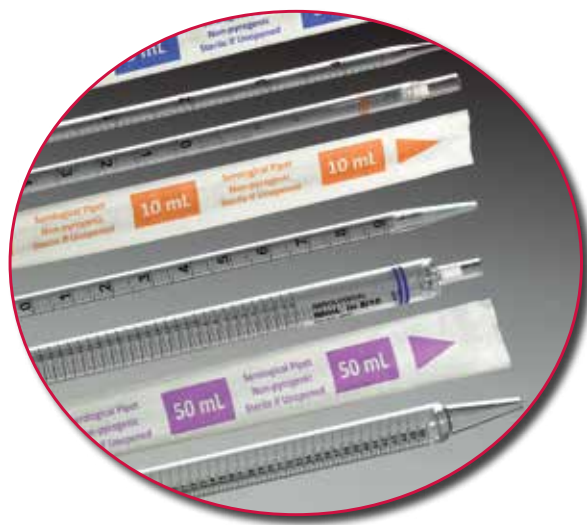


# Serological Pipets

- **Optically clear polystyrene**
- **Color coded for easy identification**
- **Individually packaged; sterile**

Thomas® brand serological pipets are manufactured from optically clear polystyrene in a class 100,000 clean room. The 1ml, 2ml and 5ml are a one piece construction and the 10ml, 25ml and 50ml are made with a durable, ultrasonically welded three piece construction. Graduations are calibrated for accurate dispensing to within  $\pm 2\%$ . Pipets are color coded for easy identification. Individually wrapped pipets are in sterile easy to open paper/plastic packaging with easy to read size indicators. The pipet plugs are a premium material which creates a filter barrier, which performs better than standard cotton plugs. Bulk packed pipets are also available in all sizes. Gamma irradiation sterilized, RNase-free, DNase-free and Non-pyrogenic.

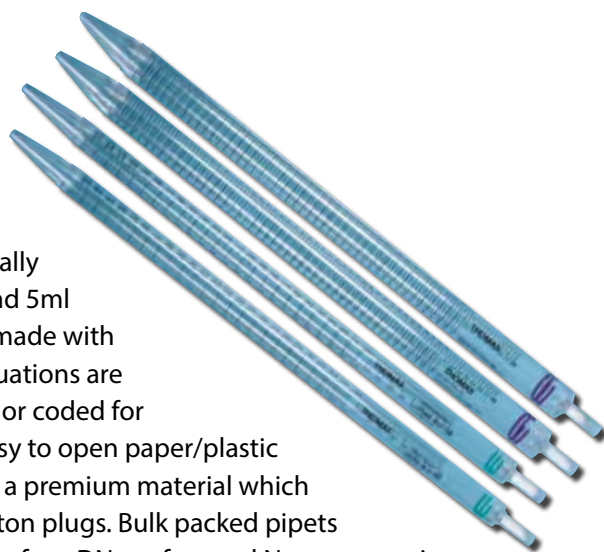

**Samples Available!**

Visit: [www.thomassci.com/samples](http://www.thomassci.com/samples)

Thomas No.	Capacity (ml)	Subdivision	Package Color Code	Case
<i>Individually Wrapped</i>				
1220Z78	1	1/100	Yellow	800
1220Z79	2	1/50	Green	600
1220Z80	5	1/10	Blue	200
1220Z81	10	1/10	Orange	200
1220Z82	25	2/10	Red	200
1220Z83	50	5/10	Purple	90
<i>Bulk</i>				
1220Z84	1	1/100	Yellow	1000
1220Z85	2	1/50	Green	800
1220Z86	5	1/10	Blue	500
1220Z87	10	1/10	Orange	500
1220Z88	25	2/10	Red	200
1220Z89	50	5/10	Purple	200